FR BE

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[www.kompernass.com](http://www.kompernass.com)

# Contents

<b>Introduction</b> .....	<b>26</b>
Information about these operating instructions. ....	26
Copyright .....	26
Intended use .....	26
Warnings .....	27
<b>Safety</b> .....	<b>27</b>
Basic safety instructions .....	27
Information on using batteries .....	29
<b>Operating components</b> .....	<b>30</b>
<b>Operation</b> .....	<b>30</b>
Package contents and transport inspection .....	30
Power supply .....	30
<b>Handling and use</b> .....	<b>31</b>
Fitting the cleaning attachment. ....	31
Cleaning the interdental spaces .....	32
<b>Cleaning</b> .....	<b>34</b>
<b>Storage/disposal</b> .....	<b>34</b>
Storage. ....	34
Disposal of the appliance. ....	34
Disposal of (rechargeable) batteries .....	35
Disposal of packaging .....	35
<b>Appendix</b> .....	<b>35</b>
Technical details .....	35
Warranty .....	35
Service .....	36
Importer .....	36

## Introduction

### Information about these operating instructions

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product.

They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please pass on all documentation incl. these operating instructions to any future owner(s) if you sell this product or give it away.

### Copyright

This documentation is protected by copyright.

Any copying or reproduction, including in the form of extracts, or any reproduction of images (even in a modified state), is permitted only with the written authorisation of the manufacturer.

### Intended use

This appliance is designed for cleaning the gaps between your teeth (interdental spaces). The appliance is not intended for any other purpose, nor for use beyond the scope described. The appliance is not intended for use in commercial or industrial environments.

The manufacturer accepts no responsibility for damage caused by failure to observe these instructions, improper use or repairs, unauthorised modifications or the use of unapproved replacement parts. The operator bears sole liability.

## Warnings

The following warning notice types are used in these operating instructions:

### CAUTION

**A warning notice at this hazard level indicates a risk of property damage.**

Failure to avoid this situation could result in material damage.

- ▶ Follow the instructions in this warning notice to prevent material damage.

### NOTE

- ▶ A note provides additional information that will assist you in using the appliance.

## Safety

This section contains important safety instructions for using the appliance.

This appliance complies with statutory safety regulations. Improper use may result in personal injury and property damage.

### Basic safety instructions

To ensure safe operation of the appliance, follow the safety guidelines set out below:

- Check the appliance for visible external damage before use. Do not operate an appliance that has been damaged or dropped.
- Cleaning attachments may be used by children and people with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been instructed how to use the appliance safely and have understood the potential risks.

- Children must not play with the appliance.
- Cleaning and user maintenance should not be performed by children unless they are under supervision.
- Do not allow children to play with packaging materials. **There is risk of suffocation.**
- The cleaning attachments ② may not be swallowed or inhaled. If this does happen, however, contact a doctor immediately.
- All repairs must be carried out by authorised specialist companies or by the Customer Service department. Improper repairs may put the user at risk. It will also invalidate any warranty claims. The hand element may not be disassembled by the user.
- Repairs to the appliance during the warranty period may only be carried out by a customer service department authorised by the manufacturer. Otherwise, no warranty claims will be held for any subsequent damages.
- Defective components must always be replaced with original replacement parts. Compliance with safety requirements can only be guaranteed if original replacement parts are used.
- Do not place any objects on the appliance.

- Never operate the appliance in the direct vicinity of open flames (e.g. candles) or in extreme environmental conditions (e.g. in an explosive atmosphere).
- If you notice any unusual noises, smoke or similar uncertain situations, remove the batteries and disconnect any attached mains adapter from the power supply and contact Customer Service (see section "Service").

## Information on using batteries

### CAUTION

- ▶ Do not throw batteries into a fire.
- ▶ Never short-circuit batteries.
- ▶ Do not attempt to recharge the batteries.
- ▶ Check the condition of the battery regularly. Leaking battery acid can cause permanent damage to the appliance.
- ▶ Special care should be taken when handling damaged or leaking batteries. Risk of acid burns! Wear protective gloves.
- ▶ Store batteries in a place that is inaccessible to children. If a battery is swallowed, seek medical attention immediately.
- ▶ If you do not intend to use the appliance for a long time, remove the battery.

## Operating components

(See fold-out page for illustrations)

- ❶ Appliance head
- ❷ Cleaning attachment
- ❸ Attachment fitting
- ❹ Hand element
- ❺ Operating button
- ❻ Battery
- ❼ Housing
- ❽ Storage box

## Operation

### Package contents and transport inspection

Check the contents of the package. The product includes the following components:

- Interdental cleaner
- 1 × 1.5 V battery
- 30 cleaning attachments
- Operating instructions

#### NOTE

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Service hotline (see section **Service**).
- ◆ Remove all parts of the appliance and the operating instructions from the carton.
- ◆ Remove all packaging material.

### Power supply

- ◆ Use a 1.5 V battery ❻ to operate the appliance. To insert the supplied battery ❻ proceed as follows:

- 1) Carefully pull the housing **7** off the hand set **4** and place it to one side.
- 2) Insert the battery **6**. Pay attention to the polarity markings given in the battery compartment.
- 3) Push the housing **7** back onto the hand set **4** until it sits firmly in place and the gap is fully closed. The anti-slip dimples on the housing **7** must be on the same side as the anti-slip dimples on the hand set **4**. The writing on the housing **7** must be on the same side as the operating button **5**.

## NOTE

- ▶ If the effectiveness of the cleaning starts to decrease, the battery **6** needs to be replaced.

## Handling and use

This section provides important information about proper handling and operation of the appliance.

### Fitting the cleaning attachment


- 1) Remove the appliance head **1** by turning it slightly so that the dot on the appliance head **1** is aligned with the  symbol on the hand set **4**. Then pull the appliance head **1** upwards and off.
- 2) Push the attachment fitting **3** into the hole of one of the cleaning attachments **2** in the storage box **8**. You will feel a slight click.
- 3) Pull the hand set **4** away from the storage box **8**. The cleaning attachment **2** comes out of the storage box **8**.
- 4) Check that the cleaning attachment **2** is firmly attached to the attachment fitting **3**. If necessary, push the cleaning attachment **2** onto the attachment fitting **3** with slightly more force.
- 5) Replacing the appliance head **1**: guide the fitted cleaning attachment **2** through the hole in the appliance head **1**. The dot on the appliance head **1** must align with the  symbol. Turn the appliance head **1** until the dot is lined up with the  symbol.

## Cleaning the interdental spaces

### NOTE

- ▶ We recommend using the appliance in front of a mirror so that you can guide the appliance more confidently and avoid injuring your gums.
- ▶ You can use this appliance every day to clean your interdental spaces. If you do so, we recommend changing the cleaning attachment **2** roughly every 3 days, as every use causes the edges of the cleaning attachment **2** to round off a little, reducing the cleaning effectiveness.

- 1) Hold the appliance vertically during the entire cleaning process.
- 2) Guide the cleaning attachment **2** from the outside into the interdental space you wish to clean. The cleaning attachment **2** is long enough: there is no need to also clean the interdental spaces from the rear.

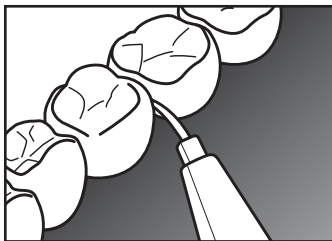- 3) Press the operating button **5**. The cleaning attachment **2** starts moving.
- 4) Move the appliance up and down to clean the entire lateral surfaces of the adjacent teeth.

### NOTE

- ▶ Do not use force when using the appliance! If you experience discomfort during one cleaning phase, continue with a different interdental space.
- 5) Once the interdental space appears sufficiently clean, release the operating button **5**. The appliance stops.

- 6) Pull the cleaning attachment **2** out of the interdental space and continue cleaning with the next interdental space.


### NOTE

- ▶ Divide your jaw into 4 quadrants. Once you have cleaned all the interdental spaces in one quadrant, rinse the cleaning attachment **2** under running water.

- 7) After completing the cleaning session, rinse the cleaning attachment **2** thoroughly under running water. Hold the appliance so that the cleaning attachment **2** is pointing downwards and the water flows over the cleaning attachment.

If you want to remove and dispose of the cleaning attachment **2**:

- 1) Remove the appliance head **1** by turning it a little so that point on the appliance head **1** is aligned with the  symbol on the hand set **4**. You can then remove the appliance head **1**.
- 2) Pull the cleaning attachment **2** off the attachment fitting **3**. If it is attached too firmly, use the removal tool on the storage box **8**.
- 3) To do this, push the cleaning attachment **2** all the way through the opening in the storage box **8** and move the hand set **4** together with the appliance head **1** in the direction of the arrow:

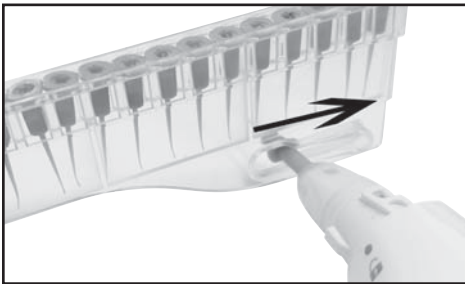This will remove the cleaning attachment **2** from the attachment fitting **3** and it will fall off.

## NOTE

- ▶ To reorder cleaning attachments **2**, please contact our service department (see section "Service"). Have the IAN number of the item to hand. You will find the IAN on the front of these operating instructions and directly under the service numbers.

## Cleaning

### CAUTION

#### Damage to the appliance!

- ▶ To avoid irreparable damage to the appliance, ensure that no moisture can penetrate the appliance during cleaning.
- Clean the surfaces of the appliance with a damp cloth. Use a normal commercial detergent and a little water to remove stubborn residue. Ensure that no water gets into the appliance.
- To clean the appliance head **1**, remove it from the hand set **4** and rinse it under running water. Allow the appliance head **1** to dry before replacing it on the hand set **4**.

## Storage/disposal

### Storage

If you plan not to use the appliance for an extended period, store it in a clean, dry place out of direct sunlight.

### Disposal of the appliance

**Never dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.**

Dispose of the appliance via an approved waste disposal company or your municipal waste disposal facility. Please comply with all applicable regulations. Please contact your waste disposal centre if you are in any doubt.

## Disposal of (rechargeable) batteries

**Used batteries must not be disposed of in the household waste.**

Consumers are legally obliged to dispose of (rechargeable) batteries at a collection point in their community/city district or at a retail store.

The purpose of this requirement is to ensure that batteries are disposed of in an environmentally friendly manner. Only dispose of batteries when they are fully discharged.

## Disposal of packaging

The packaging is made of environmentally friendly materials which you can dispose of at your local recycling centre.

## Appendix

### Technical details

Voltage supply	1.5 V battery, size AA
Protection class	IPX4

### Warranty

This appliance is provided with a 3-year warranty valid from the date of purchase. This appliance has been manufactured with care and inspected meticulously prior to delivery.

Please retain your receipt as proof of purchase. In the event of a warranty claim, please contact your Customer Service unit by telephone. This is the only way to guarantee free return of your goods.

#### NOTE

- ▶ The warranty only covers claims for material and manufacturing defects, not for damage sustained during transport, wearing parts or damage to fragile components such as switches or batteries.

This appliance is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Your statutory rights are not restricted by this warranty. The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components.

Any damage and defects present at the time of purchase must be reported immediately after unpacking.

Repairs carried out after expiry of the warranty period shall be subject to charge.

## Service

### **GB** Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: [kompernass@lidl.co.uk](mailto:kompernass@lidl.co.uk)

IAN 270480

## Importer

KOMPERNASS HANDELS GMBH

BURGSTRASSE 21

44867 BOCHUM

GERMANY

[www.kompernass.com](http://www.kompernass.com)