SILVERCREST®

www.lidl-service.com

PARAFFIN WAX BATH SPB 97 A1

GB CY

PARAFFIN WAX BATH

Operating instructions

BG

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
BG

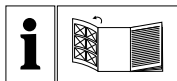GB CY

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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
B

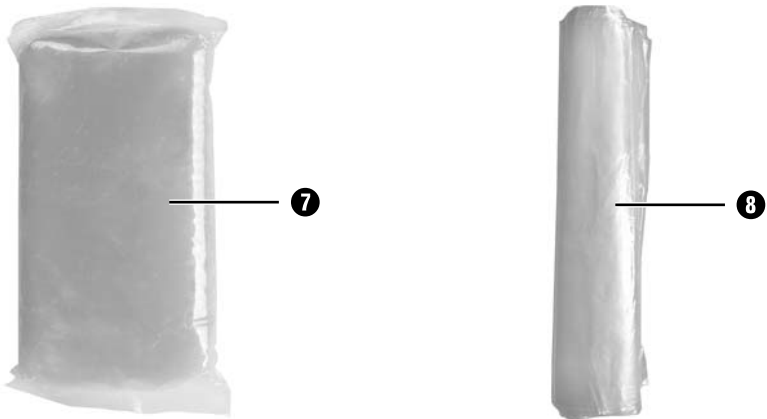Contents

Introduction	2
Intended use	2
Package contents	2
Description of the appliance	2
Technical data	2
Safety instructions	3
Why a paraffin bath?	5
Before first use	5
Preparing the appliance	5
Use of the appliance	5
Cleaning and maintenance	6
Storage	7
Disposal	7
Warranty and service	8
Importer	8

PARAFFIN WAX BATH

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product.

The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified areas of application. Please also pass on these operating instructions to any future owner.

Intended use

This appliance is intended exclusively for the heating/melting of paraffin wax. This appliance is intended solely for use in private households. Do not use it for commercial purposes.

Package contents

Paraffin wax bath

30 plastic bags

2 x 450 g paraffin wax

Operating instructions

Description of the appliance

- 1 Lid
- 2 Grille
- 3 Inner bowl
- 4 HEAT indicator lamp
- 5 Temperature control
- 6 POWER indicator lamp
- 7 Paraffin wax
- 8 Plastic bags

Technical data

Voltage: 220 - 240 V ~, 50/60 Hz

Power consumption: 97 W

Protection class: II 

Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- ▶ To avoid potential risks, damaged mains cables should be replaced by the manufacturer or the manufacturer's customer service department or a similarly qualified person.
- ▶ Use the appliance only in dry indoor areas, not outdoors.
- ▶ NEVER open the housing of the appliance.
- ▶ Do not use the appliance in the vicinity of water, such as wash basins, showers or bath tubs.
- ▶ Do not use the appliance while bathing, showering or sleeping.
- ▶ Never touch the mains plug or the mains cable with wet hands.
- ▶ Check the appliance for visible external damage before use. Do not operate an appliance that has been damaged or dropped.

WARNING! RISK OF INJURY!

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Children must not use the appliance as a plaything.
- ▶ Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.
- ▶ People who are insensitive to heat may not use this device.
- ▶ Take care when handling freshly melted paraffin wax. Risk of burns!

⚠ WARNING! RISK OF INJURY!

- ▶ Do not use the appliance if you have heart disease, skin redness, inflammation, swelling, sensitive areas, diabetes, open wounds, burns, eczema or psoriasis.
- ▶ If you are unsure whether you can safely use the appliance due to pre-existing conditions, consult a doctor prior to use.
- ▶ Do not use the appliance beneath blankets or towels.
- ▶ Do not leave the paraffin wax on your skin for longer than 20 minutes.
- ▶ Install the appliance as close as possible to an electrical power socket. Ensure that the power plug is quickly reachable in case of danger, and also that there is no danger of tripping over the cable.
- ▶ Provide a stable location for the appliance.

CAUTION! MATERIAL DAMAGE!

- ▶ Never operate this appliance when it is empty (without paraffin wax).
- ▶ Only paraffin wax may be put into the appliance. Do not add any aggressive substances or substances containing solvents!
- ▶ Always place the appliance on a heat-resistant, level surface.
- ▶ Fill the appliance no further than the MAX mark inside the appliance.
- ▶ Never heat the paraffin wax in the plastic film.
- ▶ Never heat the paraffin wax in a microwave oven, over an open flame or on a stove.
- ▶ Do not use an external timer switch or a separate remote control system to operate the appliance.
- ▶ Never use metallic, pointed or sharp objects to remove solidified paraffin wax from the appliance.

Why a paraffin bath?

Paraffin baths are used to make cracked or brittle skin smooth, supple and elastic again. They promote the blood circulation of the skin and even help to improve extremely dry skin.

The deep heat that is generated also feels very pleasant.

The heat opens the pores of the skin, allowing the active ingredients to penetrate more easily.

Before first use

- 1) Wipe the device with a damp cloth inside and out to remove any dust residue left over from transport. Dry it well so that there is no moisture left in the appliance.
- 2) Place the appliance on a level and heat-resistant surface. Also keep in mind that paraffin wax splashes may drip onto the surface. The surface should therefore be easy to clean.

Preparing the appliance

Once the appliance is set up on a suitable surface:

- 1) Remove the paraffin wax **7** from its packaging.
- 2) Put the paraffin wax **7** into the inner bowl **3**.

NOTE

- ▶ Add paraffin wax **7** to the inner bowl **3** until at least the level of the MIN marking. Do not exceed the MAX marking.

- 3) Insert the plug into a mains power socket. The POWER indicator lamp **6** illuminates.
- 4) Set the temperature control **5** to "hot". The HEAT control lamp **4** illuminates in orange to show that the appliance is heating up.
- 5) Close the lid **1**.

After about 2.5 hours the paraffin wax **7** has melted.

⚠ WARNING – RISK OF INJURY!

Never dip your hands or any other body parts into the paraffin wax **7** once the appliance has heated up to "hot". Risk of burns!

- 6) When the paraffin wax **7** has melted, remove the lid **1**.
- 7) Put the grille **2** into the inner bowl **3** and let it sink to the bottom.
- 8) Turn back the temperature control **5** until you feel a noticeable resistance when turning and the temperature control **5** is pointing to the gap on the temperature scale.
- 9) Wait until the HEAT control lamp **4** lights up green. Now the paraffin wax **7** has cooled to around 60°C and has the ideal temperature for application.

You can now start application.

Use of the appliance

NOTE

- ▶ Exterior influences such as the room temperature can affect the temperature of the paraffin wax **7**. You should therefore test to see whether the temperature is comfortable before each application. To do this, place a little paraffin wax **7** onto the inside of your wrist. If this feels pleasant, you can keep the paraffin wax at this temperature.
- ▶ If a skin forms on the paraffin wax **7**, increase the temperature slightly.
- ▶ If you wish to treat both hands at the same time you will require the assistance of a second person, since you cannot use your hands with the layer of wax on them.

You can use the paraffin bath on your hands, feet or elbows. The procedure is always the same. In this manual, the application is described using the example of the hands.

NOTE

- ▶ Depending on the condition of your hands, we recommend treatment with the paraffin bath about 1–2 times per month.

CAUTION – PROPERTY DAMAGE

- ▶ Ensure that no paraffin wax **7** gets onto your clothes during use. This can be very difficult to remove and may result in discolouration.
- 1) Remove any jewellery, such as rings or bracelets.
 - 2) Treat your hands with a rich cream, such as hand cream.
 - 3) Briefly immerse your hand in the paraffin wax **7** for about 1 - 2 seconds. Remove your hand again and wait for the paraffin wax **7** to harden.
 - 4) Repeat this process about 6 - 7 times so that a layer of paraffin wax **7** forms.

NOTE

- ▶ If you have the feeling that the temperature is too high or too low at any time, correct this using the temperature control **5**.
 - ▶ If you want to pause the heating of the appliance for a short time, turn the temperature control **5** to the left as far as the stop. The HEAT control lamp **4** goes out. This is not an adequate substitute for removing the power plug after finishing the treatment!
- 5) Insert your hand into one of the plastic bags **8**, and twist the bag **8** closed at the opening so that no air can escape.

WARNING – RISK OF INJURY!

If you pull the plastic bag **8** over your foot, do not walk around with it on. Risk of slipping.

- 6) After about 15 - 20 minutes, remove the plastic bag **8** together with the paraffin wax **7** and dispose of both in your household waste. For hygiene reasons, paraffin wax **7** that has already been applied to the body must not be reused!

NOTE

- ▶ Always remove the grille **2** from the liquid paraffin wax **7** immediately after completing the application.

NOTE

- ▶ Pure paraffin wax cannot spoil. After every use of the paraffin bath, impurities can get into the paraffin wax **7**. We therefore recommend replacing the paraffin wax **7** that has already been melted in the appliance after no more than one year.

Cleaning and maintenance

NOTE

- ▶ If the paraffin bath is only being used by one person, replace the paraffin wax **7** as soon as you notice dirt or suspended particles in the paraffin wax **7**. The appliance should also be cleaned when doing this.
- ▶ If the paraffin bath is being used by several people, the paraffin wax **7** should be changed before every treatment for hygiene reasons.
- ▶ To reorder paraffin wax **7** or plastic bags **8**, contact our Customer Service (see section "Warranty and service"). Have the IAN number of the item concerned to hand. You will find this on the cover of the operating instructions and below the service address.

To remove the paraffin wax **7** proceed as follows:

- Remove the power plug from the mains power socket.
- Leave the paraffin wax **7** to harden overnight.
- Replace the plug in the mains power socket.
- Put the lid **1** on.
- Set the temperature control **5** to "hot" and wait for approximately 10 minutes (depending on the room temperature).
- Remove the power plug from the mains power socket.

- Remove the lid ❶.
- Pull the side walls of the appliance a little apart so that the wax block detaches from the side walls.
- Press the wax block down on all sides so that you can feel whether the wax block has detached from the side walls. If the wax block does not detach, replace the power plug back in the power outlet and put the lid ❶ back on. Set the temperature control ❺ to “hot”.
- After a further 5 minutes, check whether the wax block detaches. If the wax block still does not detach, wait a little longer and check every 5 minutes whether you can detach the wax block.

NOTE

- ▶ Depending on the ambient temperature, the time that the wax requires to become soft enough to detach from the appliance can vary significantly.

Once the wax block has detached from the side walls:

- Remove the power plug from the mains power socket.
- Press on one side of the paraffin wax block. It slides upwards a little so that it can now be easily removed.
- Any liquid paraffin wax ❷ that remains in the inner bowl ❸ can be removed with a paper towel or similar.
- Then wipe the inner bowl ❸ with a cloth moistened with warm water. If necessary, add a little mild detergent onto the cloth. Afterwards, wipe the bowl with damp cloth to remove all detergent residue.
- Wipe the appliance and the lid ❶ with a moist cloth. If required, use a mild detergent on the cloth.

- Clean the grille ❷ in warm water with a mild detergent. Then rinse the grille ❷ with plenty of clean water to remove all detergent residue.
- Ensure that the appliance and all components are completely dry before using them again.

NOTE

- ▶ If any paraffin wax ❷ has dripped onto the work surface, wipe it up with a cloth moistened with warm water.

Storage

- Always clean the appliance thoroughly before storing.
- If you do not intend to use the appliance for a long period, remove all paraffin wax residue from the inner bowl ❸.
- Store the appliance in a clean and dry location with the lid ❶ on.

Disposal

You can dispose of the paraffin wax ❷ in the normal waste.

Do not dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved waste disposal company or your municipal waste disposal facility. Please comply with all applicable regulations. Please contact your waste disposal facility if you are in any doubt.

The packaging is made of environmentally friendly materials which you can dispose of at your local recycling centre.

Warranty and service

This appliance is provided with a 3-year warranty valid from the date of purchase. The appliance has been manufactured with care and inspected meticulously prior to delivery. Please retain your receipt as proof of purchase. In the event of a warranty claim, please contact your Customer Service unit by telephone. This is the only way to guarantee free return of your goods.

The warranty only covers claims for material and manufacturing defects, not for damage sustained during carriage, components subject to wear and tear or damage to fragile components such as switches or batteries. This appliance is intended solely for private use and not for commercial purposes.

The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Your statutory rights are not restricted by this warranty.

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking, but no later than two days after the date of purchase. Repairs carried out after expiry of the warranty period shall be subject to charge.

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