ELECTRIC HARD SKIN REMOVER SHE 3 B2

(GB) (IE)

ELECTRIC HARD SKIN REMOVER

Operating instructions

(FR) (BE)

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Mode d'emploi

(DE) (AT) (CH)

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(DK)

ELEKTRISK FODFIL

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(NL) (BE)

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
IAN 273036

(GB) (IE)

(DK) (BE) (NL)

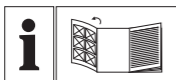GB IE

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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A**B**

Contents

Introduction	2
Information concerning these operating instructions	2
Copyright	2
Limitation of liability	2
Intended use	3
Package contents	4
Safety	5
Operating components	9
Technical data	9
General	9
Operation	10
Inserting/replacing the batteries	10
Fitting/replacing the pumice attachment	10
Use of the appliance	11
Cleaning and care	13
Storage	13
Disposal	14
Disposal of the appliance	14
Disposal of (rechargeable) batteries	14
Disposal of packaging	14
Warranty	15
Service	16
Importer	16

Introduction

Information concerning these operating instructions

Congratulations on the purchase of your new appliance. You have selected a high-quality product. The operating instructions are part of this product. They contain important information on safety, usage and disposal. Before using the product, familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified areas of application. Keep these operating instructions in a safe place. Please also pass on these operating instructions to any future owner.

Copyright

This documentation is protected by copyright.

Any copying or reproduction, including in the form of extracts, or any reproduction of images (even in a modified state), is permitted only with the written authorisation of the manufacturer.

Limitation of liability

All technical information, data and information for installation and operation contained in these operating instructions corresponds to the latest version(s) at the time of printing and, to the best of our knowledge, takes into account our previous experience and know-how.

No claims can be derived from the details, illustrations and descriptions in these instructions.

The manufacturer assumes no responsibility for damage caused by failure to observe these instructions, improper use, incompetent repairs, unauthorised modifications or for the use of non-approved replacement parts.

Intended use

This appliance is intended exclusively for the removal of hard skin (calluses). It is intended for domestic use only.

Not for commercial use!

The appliance is not intended for any other purpose, nor for use beyond the scope described.

WARNING

Danger if not used for intended purpose!

Failure to use the appliance for its intended purpose and/or any other type of use can be hazardous.

- ▶ Use the appliance only for its intended purpose.
- ▶ Observe the procedures described in these operating instructions.

No claims of any kind will be accepted for damage arising as a result of incorrect use.

The risk shall be borne solely by the user.

Package contents

The appliance is supplied with the following components as standard:

- Hand element
- Protective cap
- Coarse pumice attachment (grey)
- Fine pumice attachment (red)
- Cleaning brush
- Operating instructions
- ◆ Remove all parts of the appliance from the packaging and remove all packaging materials.

WARNING

Risk of suffocation!

- ▶ Do not allow children to play with packaging materials. There is a risk of suffocation.

NOTE

- ▶ Check the package for completeness and for signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Service hotline (see section **Service**).

Safety

⚠ WARNING – RISK OF INJURY!

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
Children must not use the appliance as a plaything. Cleaning and user maintenance tasks may not be carried out by children unless they are supervised.
- ▶ Check the appliance for visible external damage before use. Do not operate an appliance that has been damaged or dropped.
- ▶ All repairs must be carried out by authorised specialist companies or by the Customer Service department. Improper repairs can pose significant risks to the user. It will also invalidate any warranty claims.

⚠ WARNING – RISK OF INJURY!

- ▶ Extended use may cause extreme warming of the skin. Therefore, take frequent breaks. Otherwise, there is a risk of burns!
- ▶ Do not use the appliance while bathing or showering.
- ▶ Never touch the appliance with wet or moist hands.

CAUTION – PROPERTY DAMAGE

- ▶ Defective components must always be replaced with original replacement parts. Compliance with safety requirements can only be guaranteed if original replacement parts are used.
- ▶ Protect the appliance from penetration by moisture and liquids.
- ▶ Protect the appliance against knocks, dampness, dust, chemicals, extremes of temperature, etc. and keep it away from sources of heat (ovens, radiators).
- ▶ Never use the appliance under blankets, pillows, etc.

CAUTION – PROPERTY DAMAGE

- ▶ Do not store the appliance where it could fall into a bath tub or wash basin.
- ▶ NEVER submerge the appliance in water or other liquids!
- ▶ Ensure that the drive can always rotate freely. Otherwise, the appliance may overheat and be damaged beyond repair.

INFORMATION ON USING BATTERIES

- ▶ If you do not intend to use the appliance for an extended period, remove the batteries – batteries can discharge themselves and then leak.
- ▶ If the batteries leak: Put on protective gloves and clean the battery compartment with a dry cloth.
- ▶ Use **ONLY** the type of battery detailed in the Technical Data as replacements for discharged batteries.
- ▶ Do not throw the batteries into the fire. Do not short-circuit the batteries. Do not attempt to recharge the batteries. There is a risk of explosion!
- ▶ Keep batteries away from small children. Infants can put batteries into their mouths and swallow them. If a battery is swallowed, seek medical assistance **IMMEDIATELY**. There is a risk of fatal injury!

Operating components

- 1 Protective cap
- 2 Catch
- 3 Coarse pumice attachment (grey)
- 4 Release
- 5 ON/OFF switch with integrated safety button
- 6 Battery compartment
- 7 Cleaning brush
- 8 Fine pumice attachment (red)
- 9 Drive shaft
- 10 Drive shaft support

Technical data

Appliance	
Power supply	2x 1.5 V LR6/AA batteries 
Operating speeds	approx. 30 rpm
Protection class	III 

General

Hard skin (calluses) can form where the skin is exposed to particular stress due to rubbing and pressure.

If the hard skin becomes too thick, e.g. due to ill-fitting shoes or very dry skin, this can lead to pain.

If necessary, clarify the cause of excessive callus formation with your doctor.

Operation

Inserting/replacing the batteries

- 1) Open the battery compartment **6** by pressing down on the grooved part of the battery compartment cover and pushing the battery compartment cover downwards.
- 2) Remove the old batteries if inserted.
- 3) Insert new batteries as shown in the polarity diagram in the battery compartment **6**.
Always use 1.5 V AA/LR6 batteries.
- 4) Push the battery compartment cover back onto the battery compartment **6** until it clicks into place.

Fitting/replacing the pumice attachment

NOTE

- ▶ To order new pumice attachments **3/8**, contact our Service team (see section "Service"). Please ensure that you have the IAN number to hand (see cover of the operating instructions).

To remove a pumice attachment **3/8**:

- 1) Press the release **4** downwards and, at the same time, pull the side of the pumice attachment **3/8** next to the release **4** upwards.
- 2) You can then pull the pumice attachment **3/8** off the drive shaft **9**.

To fit a pumice attachment ③/⑧ (see fig. 1):

- 1) Push the pumice attachment ③/⑧ with the drive shaft support ⑩ onto the drive shaft ⑨.
- 2) Press the other side of the pumice attachment ③/⑧ downwards. Ensure that the catch ② slides correctly into the groove and clicks into place.

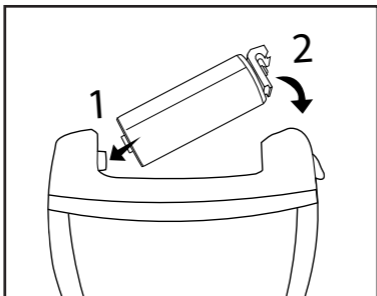
Figure 1

Use of the appliance

⚠ WARNING – RISK OF INJURY!

- ▶ Excessively long or intensive treatment can also injure deeper layers of skin. Therefore, check the treatment progress and results regularly. Diabetics in particular are less sensitive on their hands and feet.

NOTE

- ▶ The appliance will not be so effective on softened skin. Therefore, do not carry out any pretreatment in a water bath before using the appliance.
- ▶ Stop treatment IMMEDIATELY if it becomes uncomfortable or painful.

- 1) Remove the protective cap **1** from the appliance.
- 2) Select the desired pumice attachment **3/8**:
 - coarse pumice attachment **3** (grey):
for particularly thick calluses
 - fine pumice attachment **8** (red):
for mild calluses and to smooth the skin.
- 3) Fit the pumice attachment **3/8** as described in the section "Operation".
- 4) Press the integrated safety button in the On/Off switch **5** and, at the same time, push the On/Off switch **5** upwards. The pumice attachment **3/8** starts to turn.
- 5) Applying light pressure, move the appliance slowly over the hard skin that you want to remove. Move the appliance in different directions. During the treatment, do not linger in one particular place. The friction can cause unpleasant heat.
- 6) When you have finished the treatment, push the On/Off switch **5** downwards.
- 7) Remove any skin residue and dust from the pumice attachment **3/8** and the hand element using the supplied cleaning brush **7**.
- 8) Replace the protective cap **1**.

NOTE

- ▶ After the treatment, massage the skin with a skincare cream.

Cleaning and care

CAUTION – PROPERTY DAMAGE

- ▶ Do not use chemical, aggressive or abrasive cleaning agents. These can damage the surface of the appliance.
- ▶ Do not submerge the appliance or the pumice attachments **3/8** in water or other liquids.
 - ◆ Clean the appliance, especially the pumice attachments **3/8**, after every use with the cleaning brush **7**.
 - ◆ Clean the appliance and the pumice attachments **3/8** with a moist cloth. Ensure that no water or moisture can penetrate the appliance. Dry everything properly.
 - ◆ For hygienic reasons, we recommend cleaning the attachments with a clean cloth or brush moistened with a disinfectant after every use. Dry everything properly.

Storage

- ◆ Store the appliance in a dry and clean location.

Disposal

Disposal of the appliance

Do not dispose of the appliance in your normal domestic waste. This product is subject to the provisions of the European WEEE (Waste Electrical and Electronic Equipment) Directive 2012/19/EU.

Dispose of the appliance via an approved waste disposal company or your municipal waste disposal facility. Please comply with all applicable regulations. Please contact your waste disposal facility if you are in any doubt.

Disposal of (rechargeable) batteries

Used batteries may not be disposed of in the household waste. Consumers are legally obliged to dispose of (rechargeable) batteries at a collection point in their community/city district or at a retail store. The purpose of this requirement is to ensure that batteries are disposed of in an environmentally friendly manner. Only dispose of batteries when they are fully discharged.

Disposal of packaging

The packaging protects the appliance from damage during transport. The packaging materials have been selected for their environmental friendliness and ease of disposal, and are therefore recyclable.

Returning the packaging into the material cycle conserves raw materials and reduces the amount of waste generated. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

NOTE

- ▶ If possible, keep the appliance's original packaging during the warranty period so that the appliance can be packed properly for return shipment in the event of a warranty claim.

Warranty

This appliance is provided with a 3-year warranty valid from the date of purchase. The appliance has been manufactured with care and inspected meticulously prior to delivery.

Please retain your receipt as proof of purchase. In the event of a warranty claim, please contact your Customer Service unit by telephone. This is the only way to guarantee free return of your goods.

NOTE

- ▶ The warranty only covers claims for material and manufacturing defects, not for transport damage, worn parts or damage to fragile components such as buttons or rechargeable batteries.

This appliance is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Your statutory rights are not restricted by this warranty. The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components.

Any damage and defects present at the time of purchase must be reported immediately after unpacking, but no later than two days after the date of purchase.

Repairs carried out after expiry of the warranty period shall be subject to charge.

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