## ELECTRIC TOOTHBRUSH NAZ 2.4 C3

GB IE NI CY

### ELECTRIC TOOTHBRUSH

Operating instructions

DE AT CH

### ELEKTRISCHE ZAHNBÜRSTE

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GR CY

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IAN 288765

GB IE NI CY

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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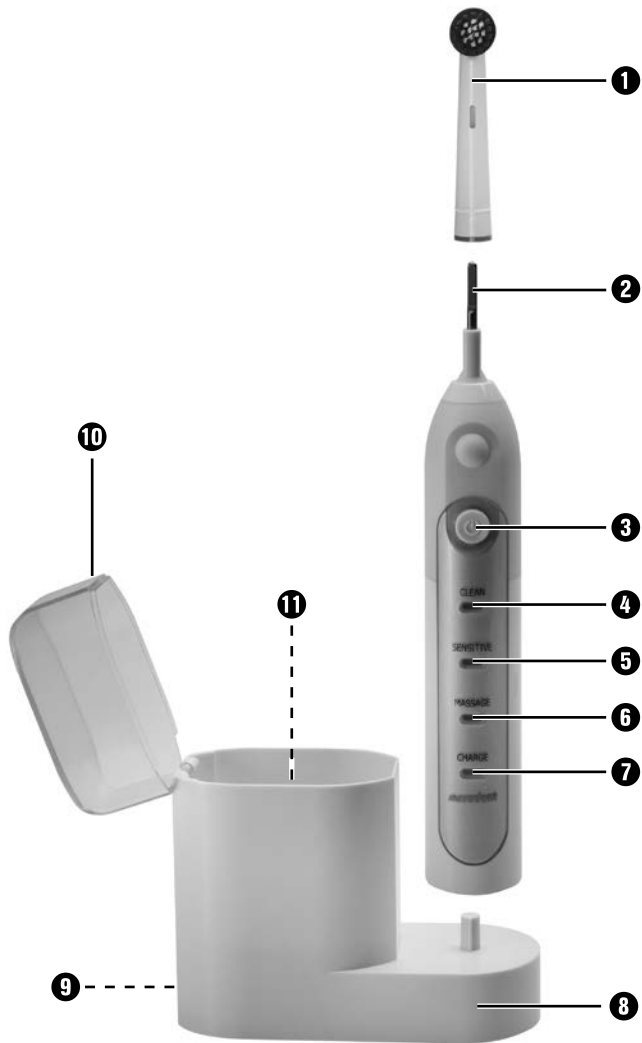
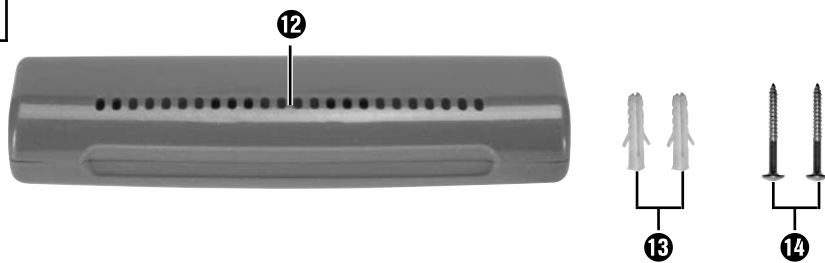
**A****B**

# Contents

<b>Introduction</b> .....	<b>2</b>
<b>Intended use</b> .....	<b>2</b>
<b>Package contents</b> .....	<b>2</b>
<b>Appliance description</b> .....	<b>3</b>
<b>Technical data</b> .....	<b>3</b>
<b>Safety information</b> .....	<b>4</b>
<b>Set-up/fitting to the wall</b> .....	<b>5</b>
<b>First use</b> .....	<b>5</b>
<b>Fitting/removing brush heads</b> .....	<b>6</b>
<b>Operation</b> .....	<b>7</b>
<b>Brushing teeth</b> .....	<b>7</b>
<b>Cleaning and maintenance</b> .....	<b>9</b>
<b>Travel lock</b> .....	<b>9</b>
<b>Travel case</b> .....	<b>10</b>
<b>Disposal</b> .....	<b>10</b>
<b>Ordering replacement parts</b> .....	<b>11</b>
<b>Kompernass Handels GmbH warranty</b> .....	<b>12</b>
Service .....	13
Importer .....	13

## Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner.

## Intended use

This appliance is intended exclusively for personal dental care at home. This appliance is intended solely for use in private households. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. Do not use it for commercial purposes.

## Package contents

- Handset
- 8 brush heads
- Travelling case
- Charger
- 2 screws
- 2 wall plugs
- Operating instructions

### **NOTE**

Check contents for completeness and damage immediately after unpacking. Contact customer service if necessary.

Remove all packaging materials and plastic films from the appliance.

## Appliance description

Figure A:

- ❶ Brush head
- ❷ Drive shaft
- ❸ On/Off switch
- ❹ "Clean" LED (orange)
- ❺ "Sensitive" LED (green)
- ❻ "Massage" LED (yellow)
- ❼ Charge indicator
- ❽ Charger
- ❾ Drill holes
- ❿ Lid
- ⓫ Brush holders

Figure B:

- ❿ Travelling case
- ⓫ Wall plugs
- ⓬ Screws

## Technical data

Power supply:	100-240 V ~ (AC), 50/60 Hz
Power consumption:	2 W
Rechargeable battery:	Li-ion battery, 600 mAh, 3.7 V
Protection type (charger):	IP X7 (protected against the ingress of water during temporary immersion)
Protection class:	II / 
Charging time:	approx. 14-16 hours (first charge)
Operating time:	approx. 60 minutes (with fully charged battery)
Protection type (handset):	IP X4 (protected against water splashes)

## Safety information

### **DANGER – RISK OF ELECTRICAL SHOCK!**

- ▶ To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
  - ▶ Ensure that no foreign objects enter the appliance.
  - ▶ In the event of malfunctions or obvious damage, please contact our customer service department.
  - ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way.  
Risk of electric shock!
  - ▶ The mains cable of this appliance cannot be replaced. If the cable is damaged, the whole appliance must be scrapped.
  - ▶ Never touch the charger if it has fallen into water. Remove the plug from the wall socket first! Risk of electric shock!
-  Do not immerse the handset or the charger in water or other liquids.

### **WARNING – RISK OF INJURY!**

- ▶ Toothbrushes may be used by children and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance should not be carried out by children, unless they are under supervision.
- ▶ Do not allow children to play with the appliance.

## CAUTION – PROPERTY DAMAGE!

- ▶ The appliance may only be recharged with the supplied charger.  
Otherwise, the appliance could be irreparably damaged!
- ▶ Operate the appliance only with the original accessories supplied.

## Set-up/fitting to the wall

### NOTE

In damp rooms, especially in the vicinity of water (bath), the charger **8** must be set up so that it cannot accidentally fall into the water.

Place the charger **8** on a clean and level surface.

Ensure the charger **8** cannot fall into water, such as a bath or sink.

You can also fit the charger permanently to the wall:

- 1) Mark two drill holes on the wall at a distance of approx. 44 mm from each other. Check to ensure that there is a mains power socket in the immediate vicinity.
- 2) Before drilling the holes, check to ensure that there are no electrical cables or water pipes behind the markings.
- 3) Use a 5 mm drill bit to drill the marked holes and push the supplied wall plugs **13** in.
- 4) Screw in the screws **14** so that the heads protrude a few millimetres. Hang the charger **8** with the drill holes **9** over the two screws **14** and ensure that the screws **14** slide into the slots. The charger **8** should lie flush against the wall. You can adjust the distance to the wall and compensate for any unevenness in the wall by screwing or unscrewing the screws.
- 5) Hang the charger **8** on the screws **14**.
- 6) Ensure that the power cable is not trapped and that it lies properly in the cut-out.

## First use

### NOTE

Before starting to use the appliance, check that it is in fault-free condition. You must not use the appliance if it has any faults.

- 1) Insert the power plug into the socket.
- 2) Before the first use, clean the new brush heads **1** with clear water.

- 3) Open the lid **10** of the brush holders **11** and place the replacement brush heads in the four slots provided. Close the lid **10** again.
- 4) Place the handset on the charger **8** so that it slides onto the pin on the base.
- 5) Allow the handset to charge up for around 14–16 hours before using it for the first time. The charging light **7** on the front of the handset lights up red during the charging process.
- 6) Use the handset repeatedly without charging the battery until the performance becomes noticeably weaker or the toothbrush stops working and the charge indicator **7** flashes red. Then recharge the battery for 14–16 hours. When the appliance is fully charged, the charge indicator **7** lights up green.
- 7) Charge and discharge the toothbrush completely 3 times. The battery will have then reached its full capacity.
- 8) You can now leave the handset in the charger **8** after every use, so as to always have full power available.

## Fitting/removing brush heads

- To remove the brush head **1**, pull it upwards and off the drive shaft **2**.
- When fitting the brush head **1**, ensure that the drive shaft **2** on the handset slides precisely into the socket on the brush head **1**.  
Push the brush head **1** down until it completely covers the drive shaft **2**.

### NOTE

You can order replacement brush heads from our Service department (see section "Ordering replacement parts").

You can also use the brush heads for the models NAZ 3.7 A1 and DAZD 3.7 LI B2 (also available from Lidl).

The toothbrush is also compatible with brush heads from the following Oral B<sup>®</sup>\* models:

Oral-B Ortho, Oral-B Precision Clean, Oral-B 3D White, Oral-B Cross Action, Oral-B Stages Power, Oral-B Sensitive, Oral-B Deep Cleaning.

\*Oral-B is a registered trademark of Procter & Gamble.

## Operation

This toothbrush has 3 different power levels:

- “Clean”: normal speed, orange “Clean” LED **4**  
For cleaning the teeth.
  - “Sensitive”: slow speed, green “Sensitive” LED **5**  
For cleaning the teeth if you have sensitive gums.
  - “Massage”: slow speed with short breaks, yellow “Massage” LED **6**  
For massaging the gums.
- Press the On/Off switch **3** once to activate the “Clean” mode. The “Clean” LED **4** lights up.
  - Press once more to switch to “Sensitive” mode. The “Sensitive” LED **5** lights up.
  - Press a third time to switch to “Massage” mode. The “Massage” LED **6** lights up.
  - To switch the toothbrush off, press the On/Off switch **3** once more.

## Brushing teeth

### **WARNING!**

Due to the unaccustomed effect of the toothbrush, there may be some initial gum bleeding. This is normal. If this continues for more than two weeks, you should consult your dentist.

- 1) Place the brush head **1** on the handset.
- 2) Moisten the brush under flowing water.
- 3) Put a little toothpaste onto the wet brush.

### **NOTE**


To avoid spray, do not switch on the appliance until it is in your mouth.

- 4) Use the On/Off switch **3** to select the desired operating mode.

### **NOTE**

This toothbrush has a 30-second interval timer, so that all four jaw quadrants can be cleaned evenly.

- 5) Divide your mouth into four parts (quadrants): upper left, upper right, lower left and lower right.

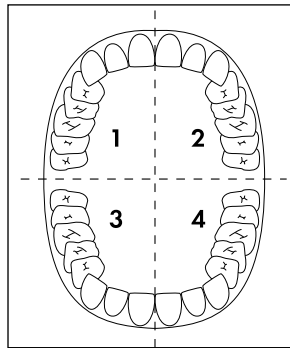- 6) Begin in one quadrant:  
Start with the brush head **1** on your gums, and move the brush upwards/ downwards from the gums to the crown of the tooth.  
Continue in this manner on both the inside and the outside of the tooth. Start again on the gums for each individual tooth.  
Do not forget the back surfaces of the rear molars.  
To clean the chewing surfaces, guide the brush over the chewing surfaces with a little pressure.
- 7) After 30 seconds, the toothbrush stops briefly to signal that you should start cleaning the next quadrant. Clean all jaw quadrants using the same technique.
- 8) After two minutes, the toothbrush stops briefly 3 times. This means that you have complied with the dental recommendation of two minutes cleaning.
- 9) Press the On/Off switch **3** repeatedly until the toothbrush stops.
- 10) Then take the toothbrush out of your mouth.
- 11) After use, clean the brush head **1** thoroughly under running water.
- 12) Remove the brush head **1** from the handset.  
To clean the drive shaft **2** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **2** is pointing downwards at an angle.

#### NOTE

Massage your gums regularly. This promotes blood circulation and helps to keep them healthy. Set the mode to "Massage" and guide the moistened brush over your gums with a light pressure.

## Cleaning and maintenance

### DANGER – RISK OF ELECTRIC SHOCK!

- ▶ Always remove the plug from the charger **8** before cleaning the appliance!

-  Do not immerse the handset or the charger **8** in water or other liquids!

### CAUTION – PROPERTY DAMAGE!

- ▶ Do not use detergents or solvents. They can damage the plastic surfaces.

### NOTE

NEVER open the housing of the charger **8** or the toothbrush! They cannot be resealed in a watertight manner.

- 1) Unplug the charger **8** from the power socket.
- 2) Clean the charger **8** inside and out with a soft, lightly moistened cloth.
- 3) Remove the brush head **1** from the handset. After use, clean the brush heads **1** thoroughly under running water.  
We recommend cleaning the brush heads **1** in a standard commercial denture cleaning fluid every 2–4 weeks. A new brush head **1** should be fitted after about 3 months, and the brush head must always be replaced when the bristles start to bend outwards.
- 4) To clean the drive shaft **2** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **2** is pointing downwards at an angle.

## Travel lock

This toothbrush is equipped with a travel lock . You can lock the On/Off button **3** to prevent accidental starting of the toothbrush:

Press and hold the On/Off button **3** until the charge indicator **7** flashes green 5 times.

The On/Off button **3** is now locked.

To unlock the switch, press and hold the On/Off button **3** until the charge indicator **7** lights up once for a few seconds.

### NOTE

- ▶ If you try to switch on the toothbrush while the travel lock is activated, the charge indicator **7** will flash 5 times. The toothbrush will not start.

## Travel case

If you wish to take the toothbrush with you on a trip, you can pack both the hand element and two brush heads ❶ in the travelling case ❷. When the battery is fully charged, you can use the toothbrush for 60 minutes without recharging. This is a guideline value only. Depending on the pressure used when brushing your teeth, the maximum operating time may vary.

## Disposal

**Never dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.**

Dispose of the appliance through an approved disposal facility or your municipal waste facility. Please comply with all applicable regulations. Please contact your waste disposal centre if you are in any doubt.

### NOTE

- ▶ The integrated rechargeable battery in this appliance cannot be removed for disposal.

The packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (b) and numbers (a) with the following meanings:


1–7: Plastics

20–22: Paper and cardboard


80–98: Composites

## Ordering replacement parts

You can order the following replacement parts for the product NAZ 2.4 C3:

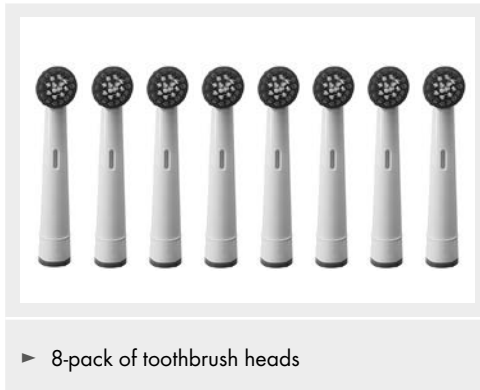Order the replacement parts via the Service hotline (see section “Service”) or simply visit our website at [www.kompernass.com](http://www.kompernass.com).
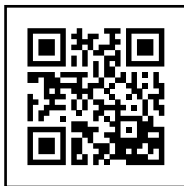
### NOTE

- Have the IAN number of the appliance ready (can be found on the cover of these operating instructions) when you place your order.

## Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

### Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

### Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

### Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

## Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on [www.lidl-service.com](http://www.lidl-service.com).

## Service

### **GB** Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: [kompernass@lidl.co.uk](mailto:kompernass@lidl.co.uk)

### **IE** Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: [kompernass@lidl.ie](mailto:kompernass@lidl.ie)

### **CY** Service Cyprus

Tel.: 8009 4409

E-Mail: [kompernass@lidl.com.cy](mailto:kompernass@lidl.com.cy)

IAN 288765

## Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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