[www.lidl-service.com](http://www.lidl-service.com)

## ZAHNPOLIERER-SET DENTAL POLISHING KIT SET DE POLISSAGE DENTAIRE NZPS 3 A I

(DE) (AT) (CH)

### ZAHNPOLIERER-SET

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(FR) (BE)

### SET DE POLISSAGE DENTAIRE

Mode d'emploi

(CZ)

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(PT)

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(GB) (IE)

### DENTAL POLISHING KIT

Operating instructions

(NL) (BE)

### TANDENPOLIJSTSET

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
(ES)

### SET DE PULIDO DENTAL

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IAN 291927

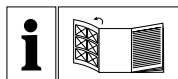DE AT CH

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# Contents

<b>Introduction</b> .....	<b>14</b>
Information about these operating instructions .....	14
Copyright .....	14
Intended use .....	14
Warnings .....	14
<b>Safety</b> .....	<b>14</b>
Basic safety instructions .....	15
Information on using batteries .....	16
<b>Operating elements</b> .....	<b>17</b>
<b>Operation</b> .....	<b>17</b>
Package contents and transport inspection .....	17
Power supply .....	18
<b>Handling and use</b> .....	<b>18</b>
Attaching the polishing cap .....	18
Removing the polishing cap .....	18
Polishing teeth .....	18
<b>Using the accessories</b> .....	<b>19</b>
Dental mirror .....	19
Interdental cleaner .....	19
Plaque and stain remover .....	19
<b>Cleaning</b> .....	<b>20</b>
<b>Storage/disposal</b> .....	<b>20</b>
Storage .....	20
Disposal of the appliance .....	20
Disposal of (rechargeable) batteries .....	20
Disposal of packaging .....	20
<b>Appendix</b> .....	<b>21</b>
Technical details .....	21
Kompermass Handels GmbH warranty .....	21
Service .....	23
Importer .....	23
Ordering replacement parts .....	24

## Introduction

### Information about these operating instructions

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please pass on all documentation incl. these operating instructions to any future owner(s) if you sell this product or give it away.

### Copyright

This documentation is protected by copyright.

Any copying or reproduction, including extracts thereof, and the reproduction of images (even in a modified state), is only permitted with the written authorisation of the manufacturer.

### Intended use

This appliance is designed for the removal of stains and discolouration from the surface of the teeth. The appliance is not intended for use in commercial or industrial environments.

The manufacturer accepts no responsibility for damage caused by failure to observe these instructions, improper use or repairs, unauthorised modifications or the use of unapproved replacement parts. The operator bears sole liability.

### Warnings

The following warnings are used in these operating instructions:

#### CAUTION

**A warning notice of this hazard level indicates a risk of property damage.**

Failure to avoid this situation could result in property damage.

- ▶ Follow the instructions in this warning notice to prevent material damage.

#### NOTE

- ▶ A note provides additional information that makes handling the appliance easier for you.

## Safety

This section contains important safety instructions for handling the appliance.

This appliance complies with statutory safety regulations. Improper use may result in personal injury and property damage.

## Basic safety instructions

To ensure safe operation of the appliance, follow the safety guidelines set out below:

- Check the appliance for visible external damage before use. Do not use an appliance that has been damaged or dropped.
- Cleaning and user maintenance should not be performed by children unless they are under supervision.
- Tooth polishers may be used by children and people with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been instructed how to use the appliance safely and have understood the potential risks.
- Do not allow children to use the appliance as a toy.
- Do not allow children to play with packaging materials.

### **There is a risk of suffocation.**

- The polishing caps **2** and the interdental cleaner attachment **13** must not be swallowed or inhaled. If this does happen, however, contact a doctor immediately.
- All repairs must be carried out by authorised specialist companies or by the Customer Service department. Improper repairs may put the user at risk. It will also invalidate any warranty claims. The hand element may not be disassembled by the user.
- Repairs to the appliance during the warranty period may only be carried out by a customer service department authorised by the manufacturer. Otherwise, no warranty claims will be held for any subsequent damages.
- Defective components must always be replaced with original replacement parts. Compliance with safety requirements can only be guaranteed if original replacement parts are used.
- Do not place any objects on the appliance.

- Never operate the appliance in the direct vicinity of open flames (e.g. candles) or in extreme environmental conditions (e.g. in an explosive atmosphere).
- If you notice any unusual noises, smoke or similar unusual situations, remove the batteries and disconnect any attached mains adapter from the power supply and contact Customer Service (see section "Service").
- This appliance is not suitable for children under the age of 3 years as it contains small parts and there is a danger of swallowing.
- The appliance and accessories may only be used by one person.
- Do not use the appliance on milk teeth.

### **Information on using batteries**

#### **CAUTION**

- ▶ Do not throw batteries into a fire.
- ▶ Never short-circuit batteries.
- ▶ Do not attempt to recharge the batteries.
- ▶ Check the condition of the battery regularly. Leaking battery acid can cause permanent damage to the appliance.

#### **CAUTION**

- ▶ Store batteries out of reach of children. If a battery is swallowed, seek medical attention immediately.
- ▶ If you do not intend to use the appliance for a long time, remove the battery.
- ▶ Special care should be taken when handling damaged or leaking batteries. Risk of acid burns! Wear protective gloves.

## Operating elements

(See fold-out page for illustrations)

Figure A:

- ❶ Attachment thread
- ❷ Polishing cap
- ❸ Operating button 
- ❹ Housing
- ❺ Hand element
- ❻ Catch 

Figure B:

- ❼ Dental mirror
- ❽ Interdental cleaner
- ❾ Plaque remover
- ❿ Stain remover
- ⓫ Batteries
- ⓬ Travelling case
- ⓭ Interdental cleaner attachment

## Operation

### Package contents and transport inspection

Check the contents of the package.

The product includes the following components:

- Tooth polisher
- 2 × 1.5 V batteries
- 4 polishing caps
- Dental mirror
- Interdental cleaner
- Plaque remover
- Stain remover
- Operating instructions

#### NOTE

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Service hotline (see section **Service**).

- ◆ Remove all parts of the appliance and the operating instructions from the carton.
- ◆ Remove all packaging material.

## Power supply

- ◆ Use two 1.5 V batteries ❶ to operate the appliance. To insert the supplied batteries ❶ proceed as follows:
  - 1) Carefully twist the housing ❷ anticlockwise off the hand element ❸ and place it to one side.
  - 2) Insert the batteries ❶. Pay heed to the correct polarity markings given on the housing ❷.
  - 3) Push the housing ❷ back onto the hand element ❸. Turn the housing ❷ clockwise to tighten it.

### NOTE

- ▶ If the polishing effect diminishes, the batteries ❶ must be replaced.

## Handling and use

This section provides important information about proper handling and operation of the appliance.

### Attaching the polishing cap

- ◆ Place the shaft of the polishing cap ❷ onto the attachment thread ❶. Turn the polishing cap ❷ clockwise while pressing down the catch  ❸ to tighten it.

### Removing the polishing cap

#### NOTE

- ▶ Remove the polishing cap ❷ when it no longer delivers the required polishing result.
- ◆ Press down the catch  ❸ while turning the polishing cap ❷ anticlockwise.

## Polishing teeth

#### NOTE

- ▶ We recommend using the appliance in front of a mirror so that you can guide the appliance more confidently and avoid injuring your gums.
- ▶ It is normally sufficient to use the appliance twice a week. Use of the appliance does not replace daily brushing with a toothbrush.
- ▶ Do not use appliance on your gums.
  - 1) Add a little standard commercial toothpaste or a special polishing toothpaste to the polishing cap ❷. Use the polishing cap ❷ to spread the toothpaste onto the area to be polished.
  - 2) Press the operating button  ❸. The polishing cap ❷ starts moving.
  - 3) Move the appliance in a circular motion up and down to polish the tooth. Polish each tooth for approx. 20–30 seconds.
  - 4) Release the operating button  ❸ when you want to stop cleaning. The appliance stops.

## NOTE

- ▶ Never apply force when using the appliance. If you experience discomfort during one polishing phase, move on to a different area.

## Using the accessories

## NOTE

- ▶ Do not use the accessories on your gums.

### Dental mirror

Use the dental mirror **7** to discover stains or plaque in hard-to-reach areas. Preferably, use the dental mirror **7** in combination with the appliance and the other accessories.

### Interdental cleaner

Use the interdental cleaner **8** to remove contamination from between the teeth.

- 1) Guide the tip of the interdental cleaner **8** into the interdental space you want to clean from the outside. As the tip is long enough, you do not need to clean the interdental spaces from the behind as well.
- 2) Move the interdental cleaner **8** back and forth to clean the entire side surface of the neighbouring teeth.

### Plaque and stain remover

Use the plaque remover **9** and the stain remover **10** to remove stains and plaque from the teeth.

- 1) Remove the protective cap.
- 2) Guide the tip of the plaque remover **9** or the stain remover **10** onto the affected area. Start with the part of the tooth closest to your gums. Now guide the tip towards to end of the tooth. Repeat the process until the plaque or stain has been successfully removed.
- 3) Replace the protective cap.

## Cleaning

### CAUTION

#### Damage to the appliance!

- ▶ To avoid irreparable damage, ensure that no moisture penetrates the appliance during cleaning.
- Clean the surfaces of the appliance and the accessories with a damp cloth moistened with a standard commercial denture cleaning fluid.
- To clean the polishing cap ②, remove it as described in the section "Removing the polishing cap". Rinse the polishing cap ② under running water. Allow the polishing cap ② to dry before screwing it back onto the attachment thread ①. We recommend cleaning the polishing caps ② in a standard commercial denture cleaning fluid every 2–4 weeks.

## Storage/disposal

### Storage

If you plan not to use the appliance for a long period, store it in a clean, dry place away from direct sunlight, preferably in the travelling case ⑫.

### Disposal of the appliance

**Never dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.**

Dispose of the appliance through an approved disposal facility or your municipal waste facility. Please comply with all applicable regulations. Please contact your waste disposal centre if you are in any doubt.

### Disposal of (rechargeable) batteries

**Used batteries must not be disposed of in the household waste.**

Consumers are legally obliged to dispose of (rechargeable) batteries at a collection point in their community/city district or at a retail store. The purpose of this requirement is to ensure that batteries are disposed of in an environmentally friendly manner. Only dispose of batteries when they are fully discharged.

### Disposal of packaging

The packaging is made of environmentally friendly materials which you can dispose of at your local recycling centre.

## Appendix

### Technical details

Power supply	2 x 1.5 V batteries type AA, LR6, DC ---
Protection type	IPX4

### Kompnass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

#### Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

#### Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

#### Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

### **Warranty claim procedure**

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on [www.lidl-service.com](http://www.lidl-service.com).

## Service

### GB Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: [kompernass@lidl.co.uk](mailto:kompernass@lidl.co.uk)

### IE Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: [kompernass@lidl.ie](mailto:kompernass@lidl.ie)

IAN 291927

## Importer

Please note that the following address is not the service address.

Please use the service address provided in the operating instructions.

KOMPERNASS HANDELS GMBH

BURGSTRASSE 21

DE-44867 BOCHUM

GERMANY

[www.kompernass.com](http://www.kompernass.com)

## Ordering replacement parts

You can order the following replacement parts for product NZPS 3 A1:

- ▶ Set: 5 polishing caps

- ▶ Set: 5 interdental cleaner attachments

Order the replacement parts via the Service Hotline (see section "Service") or simply visit our website at [www.kompernass.com](http://www.kompernass.com).

### NOTE

- ▶ Have the IAN number of the appliance ready (can be found on the cover of these operating instructions) when you place your order.

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