www.lidl-service.com

SONIC TOOTHBRUSH NSLP 3.7 A1

GB IE NI CY

SONIC TOOTHBRUSH

Operating instructions

DE AT CH

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GR CY

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GB IE NI
CY

GB IE NI CY

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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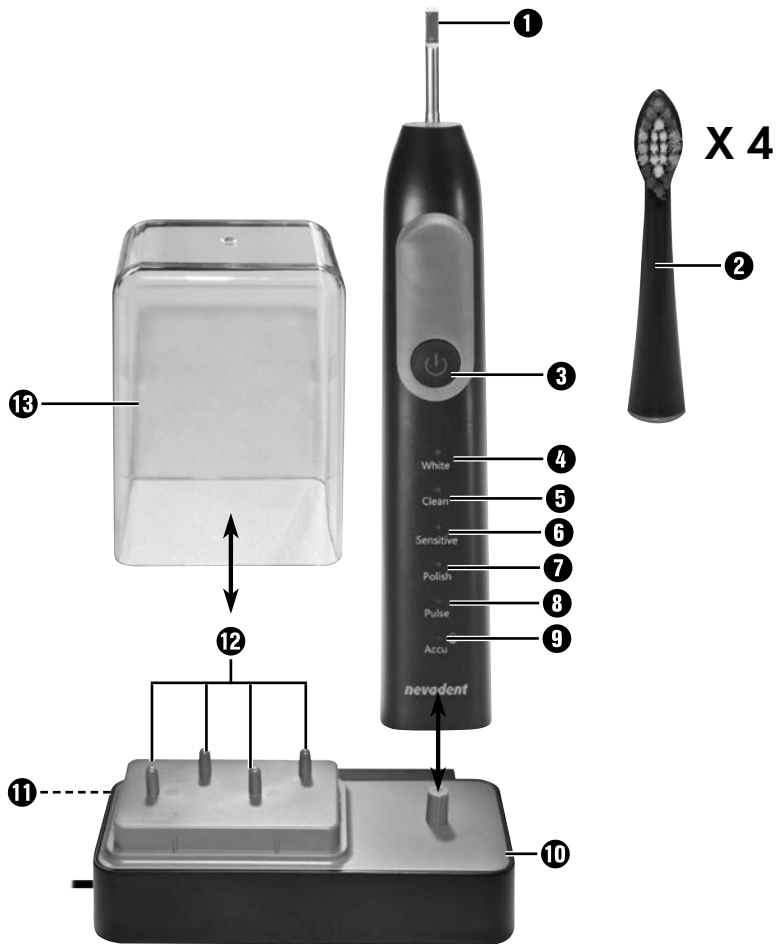
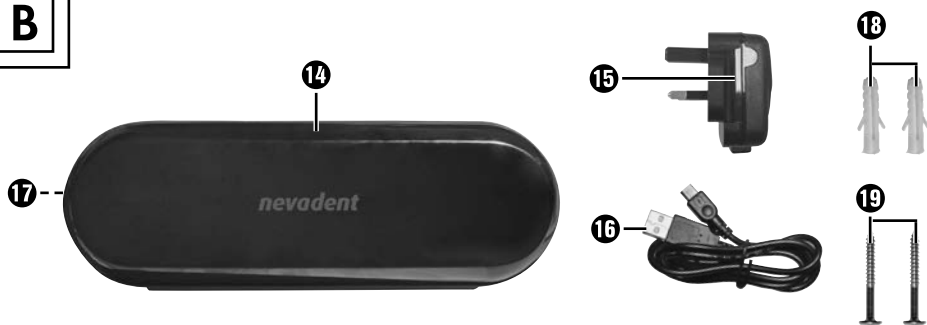
A**B**

Contents

Introduction	2
Intended use	2
Package contents	2
Appliance description	3
Technical details	4
Safety information	5
Before first use	6
Set-up/wall bracket	6
First use	7
Fitting/removing brush heads	7
Use	8
Brushing teeth	9
Cleaning	11
Travel case	12
Travel lock	12
Disposal	13
Disposal of the appliance	13
Disposal of packaging	13
Mains adapter disposal	13
Ordering replacement parts	14
Kompernass Handels GmbH warranty	15
Service	16
Importer	16

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Keep these operating instructions as a reference and store them near the product. Please also pass these operating instructions on to any future owner.

Intended use

This appliance is intended exclusively for personal dental care at home.

This appliance is intended solely for use in private households. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. Do not use it for commercial purposes.

Package contents

- Handset
- 4 brush heads
- Charging station
- USB charging cable
- Mains adapter
- 2 screws
- 2 wall plugs
- Travel case
- Operating instructions

Appliance description

Figure A:

- ① Drive shaft
- ② Brush head
- ③ On/Off button
- ④ "White" LED
- ⑤ "Clean" LED
- ⑥ "Sensitive" LED
- ⑦ "Polish" LED
- ⑧ "Pulse" LED
- ⑨ "Accu" LED
- ⑩ Charger
- ⑪ Screw holes
- ⑫ Replace brush head holders
- ⑬ Lid

Figure B:

- ⑭ Travel case
- ⑮ Mains adapter
- ⑯ USB charging cable
- ⑰ Charging socket
- ⑱ Wall plugs
- ⑲ Screws

Technical details

Mains adapter 15	
Type	GQ06-050120-AB
Input voltage	100–240 V ~ (AC current), 50/60 Hz
Power consumption	max. 0.3 A
Output voltage	5 V \equiv (DC)
Output current	1.2 A
Protection class	II / 
Short-circuit proof safety transformer	
Switched-mode power supply	
Polarity USB	P1   P4 
Charger 10	
Input voltage	100–240 V ~, 50/60 Hz
Power consumption	3.5 W
Protection class	II / 
Protection type	IPX7 Protection against temporary submersion
Handset	
Battery	3.7 V/1200 mAh Li-ion
Protection type	IPX7 Protection against temporary submersion
Travel case 14	
Input voltage	5 V \equiv
Output current	1200 mA
Protection class	

Safety information

WARNING! RISK OF ELECTRIC SHOCK!

- ▶ To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
- ▶ Ensure that no foreign objects enter the appliance.
- ▶ In the event of malfunctions or obvious damage, please contact our customer service department.
- ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way.
Risk of electric shock!
- ▶ The mains cable on this appliance cannot be replaced. If the mains cable is damaged, the whole appliance must be scrapped.
- ▶ Never touch the charger if it has fallen into water. Remove the plug from the wall socket first! Risk of electric shock!

 Do not submerge the handset or the charger in water or other liquids.

WARNING! RISK OF INJURY!

- ▶ Toothbrushes may be used by children and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance should not be carried out by children, unless they are under supervision.
- ▶ Do not allow children to use the appliance as a toy.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not use the "Polish" mode on milk teeth.

CAUTION! PROPERTY DAMAGE!

- ▶ The electric toothbrush may only be recharged with the supplied charger.
Otherwise, the appliance could be irreparably damaged!
- ▶ Operate the appliance only with the original accessories supplied.

Only the use mains adapter/charger/appliance indoors.

Before first use

- 1) Clean the appliance as described in the section "Cleaning" to remove any transport dust. Dry it thoroughly.
- 2) Place the appliance on a stable and level surface.

Set-up/wall bracket

NOTE

- ▶ The charger **10** may only be operated free-standing if used in dry rooms.
- ▶ In damp rooms, especially in the vicinity of water (bath), the charger **10** must be firmly fixed to the wall to protect it from falling into water.

- 1) Use the screw holes **11** of the charger **10** to mark the spacing for the drill holes on the wall. Check to ensure that there is a mains power socket in the immediate vicinity.
- 2) Before drilling the holes, check to ensure that there are no electrical cables or water pipes behind the markings.
- 3) Use a 5 mm drill bit to drill the marked holes and push the wall plugs **18** in.
- 4) Affix the wall bracket **10** using the screws **19** provided.
- 5) Ensure that the power cable is not trapped and that it lies properly in the cut-outs.

First use

NOTE

Before starting to use the appliance, check that it is in fault-free condition. If not, it must not be used.

- 1) Insert the plug of the charger **10** into a mains power socket.
- 2) Before first use, clean the new brush heads **2** with fresh water.
- 3) Put the replacement brush heads into the replacement brush head holders **12**. Then put the lid **13** on.
- 4) Place the handset on the charger **10** so that it slides onto the pin on the base.
- 5) Allow the handset to charge for around 18 hours before using it for the first time. The "Accu" LED **9** on the front of the handset lights up orange when the battery is under 50% charge capacity and green at more than 50% capacity during the charging process. Once the charging process is complete, the "Accu" LED **9** lights up green permanently.
- 6) Use the handset repeatedly without recharging the battery until the toothbrush pulses three times or stops working and the "Accu" LED **9** flashes red. Then, recharge the battery for 18 hours. Once the appliance is fully charged, the "Accu" LED **9** lights up green permanently.
- 7) Completely charge and discharge the toothbrush 3 times. The battery will have then reached its full capacity.
- 8) You can now leave the handset in the charger **10** after every use, so as to always have full power available.

Fitting/removing brush heads

- To remove the brush head **2**, pull it upwards and off the drive shaft **1**.
- When fitting the brush head **2**, ensure that the drive shaft **1** on the handset slides precisely into the socket on the brush head **2**. Push the brush head **2** down until it completely covers the drive shaft **1**.

Use

This toothbrush has 5 different power levels:

- "White": high speed, for especially thorough cleaning of the teeth; the "White" LED **4** lights up
- "Clean": normal speed, for daily cleaning, "Clean" LED **5** lights up.
- "Sensitive": slow speed, for people with sensitive teeth and gums, the "Sensitive" LED **6** lights up
- "Polish": highest speed, for polishing teeth, the "Polish" LED **7** lights up

NOTICE ON POLISHING

- ▶ To polish the teeth, add a little dental polishing paste to the brush head **2**. Apply the brush head **2** directly to the tooth. When polishing your teeth, do **not** start at the gums like during normal brushing. Move the appliance in a circular motion up and down to polish the tooth. Polish each tooth for approx. 20–30 seconds.

- ▶ It is normally sufficient to polish your teeth twice a week.

- "Pulse": normal speed with short breaks, to massage the gums, the "Pulse" LED **8** lights up.
- Press the On/Off button **3** once to activate the "White" mode. The "White" LED **4** and the "Accu" LED **9** light up.
- Press the On/Off button **3** once to switch up through various modes in the order described. The respective operating LED and the "Accu" LED **9** light up.
- Press the On/Off button **3** again while the toothbrush is in "Pulse" mode to switch off the toothbrush.

NOTE

If you operate one of the modes for longer than 10 seconds and then press the On/Off button **3**, the toothbrush will switch off.

Next time you switch it on, the toothbrush will automatically start in the mode you were using last.

To change the mode again, press the On/Off button **3** repeatedly until the appliance switches to the required mode as described above.

Brushing teeth

⚠ WARNING!

Due to the unaccustomed effect of the toothbrush, there may be some initial gum bleeding. This is normal. If this continues for more than two weeks, you should consult your dentist.

- 1) Place the brush head ② on the handset.
- 2) Moisten the brush under flowing water.
- 3) Put a little toothpaste onto the wet brush.

NOTE


To avoid spray, do not switch on the appliance until it is in your mouth.

- 4) Use the On/Off button ③ to select the desired operating mode.

NOTE

This toothbrush has a 30-second interval timer, so that all four jaw quadrants can be cleaned evenly.

- 5) Divide your mouth into four parts (quadrants): upper left, upper right, lower left and lower right.

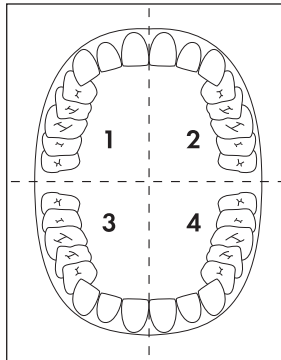- 6) Begin in one quadrant:
Place the brush head ② on your gums, and move the brush upwards/ downwards from the gums to the crown of the tooth.
Continue in this manner on both the inside and the outside of the tooth.
Start again on the gums for each individual tooth.
Do not forget the back surfaces of the rear molars.
To clean the chewing surfaces, guide the brush over the chewing surfaces with a little pressure.

- 7) After 30 seconds, the toothbrush stops briefly to signal that you should start cleaning the next quadrant. Clean all jaw quadrants using the same technique.
- 8) After two minutes the toothbrush stops. This means that you have complied with the recommendation by dental professionals of two minutes' cleaning.
- 9) Now, you can take the toothbrush out of your mouth.
- 10) After use, clean the brush head **2** thoroughly under running water.
- 11) Remove the brush head **2** from the handset.
To clean the drive shaft **1** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **1** is pointing downwards at an angle.

NOTE

Massage your gums regularly. This promotes blood circulation and helps to keep them healthy. Set the mode to "Pulse" and guide the moistened brush over your gums while applying a light pressure.

Cleaning

DANGER – ELECTRIC SHOCK!

- ▶ Always remove the plug **10** from the power socket before cleaning the appliance! There is a risk of electric shock!
-  Never immerse the handset or the charger **10** in water or other liquids!

WARNING – PROPERTY DAMAGE!

- ▶ Do not use detergents or solvents. They can damage the plastic surfaces.

NOTE

Never open the housing of the charger **10** or the toothbrush!
They cannot be resealed in a watertight manner.

- 1) Clean the charger **10** only with a soft, lightly moistened cloth.
- 2) Remove the brush head **2** from the handset. After use, clean the brush heads **2** thoroughly under running water.
We recommend cleaning the brush heads **2** in a standard commercial denture cleaning fluid every 2–4 weeks.
The brush heads **2** are equipped with so-called indicator bristles (the outer, light-blue bristles). These bristles lose their colour and thus show that the brush head **2** needs to be changed. Change the brush head **2** when the indicator bristles lose their colour.
- 3) To clean the drive shaft **1** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **1** is pointing downwards.
- 4) Clean the travel case **14** with a soft and slightly damp cloth.

Travel case

If you wish to take the toothbrush with you on a trip, you can pack the hand element and two brush heads **2** in the travel case **14**. You can also charge the toothbrush using the USB charging cable **16**:

- 1) Place the toothbrush in the appropriate slot in the travel case **14**.
- 2) Connect the USB charging cable **16** to the appropriate charging socket **17** on the travel case **14**.
- 3) Connect the other end of the USB charging cable **16** to the mains adapter **15**.
- 4) Connect the mains adapter **15** to a power socket. The flashing "Accu" LED **9** indicates that the charging process has started.

NOTE

You can use any standard commercial mains adapter (USB charger) with a suitable USB charging cable to charge the appliance. Note the technical data provided in these operating instructions.

Travel lock

This toothbrush is equipped with a travel lock. You can lock the On/Off button **3** to prevent accidental starting of the toothbrush:

Press and hold the On/Off button **3** until the "Accu" LED **9** flashes green 5 times. The On/Off button **3** is now locked.

To unlock the button, press and hold the On/Off button **3** until the "Accu" LED **9** lights up once for a few seconds.

NOTE

If you try to switch on the toothbrush while the travel lock is activated, the "Accu" LED **9** will flash 5 times. The toothbrush will not start.

Disposal

Disposal of the appliance

Never dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance through an approved disposal facility or your municipal waste facility. Please comply with all currently applicable regulations. Please contact your local waste disposal centre if you are in any doubt.

NOTE

- ▶ The integrated rechargeable battery in this appliance cannot be removed for disposal.

Disposal of packaging

The packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (b) and numbers (a) with the following meanings:

1-7: Plastics

20-22: Paper and cardboard

80-98: Composites

Mains adapter disposal

Do not dispose of the mains adapter in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the mains adapter via an approved disposal centre or at your community waste facility. Please comply with all currently applicable regulations. Please contact your local waste disposal centre if you are in any doubt.

Ordering replacement parts

You can order the following replacement parts for this product:

- ▶ 8-pack of toothbrush heads (white)

- ▶ 8-pack of toothbrush heads (black)

Order the replacement parts via the service hotline (see section "Service") or simply visit our website at www.kompernass.com.

NOTE

- ▶ Have the IAN number of the appliance ready (can be found on the cover of these operating instructions) when you place your order.

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on www.lidl-service.com.

Service

GB Service Great Britain
Tel.: 0871 5000 720 (£ 0.10/Min.)
E-Mail: kompennass@lidl.co.uk

IE Service Ireland
Tel.: 1890 930 034
(0,08 EUR/Min., (peak))
(0,06 EUR/Min., (off peak))
E-Mail: kompennass@lidl.ie

CY Service Cyprus
Tel.: 8009 4409
E-Mail: kompennass@lidl.com.cy

IAN 291945

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

KOMPENASS HANDELS GMBH
BURGSTRASSE 21
DE-44867 BOCHUM
GERMANY
www.kompennass.com

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