[www.lidl-service.com](http://www.lidl-service.com)

## SONIC TOOTHBRUSH NSZB 3.7 B3

GB IE NI

### SONIC TOOTHBRUSH

Operating instructions


DE AT CH

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IAN 303047

GB IE NI


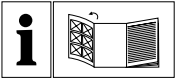
GB IE NI

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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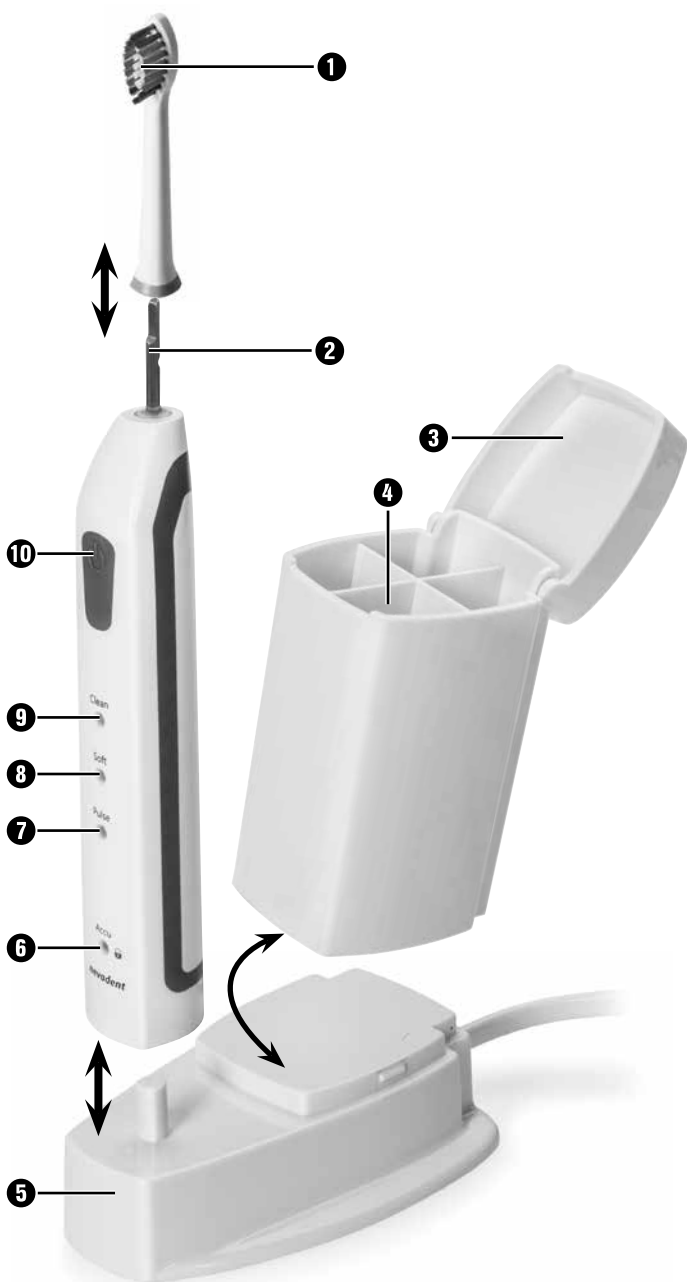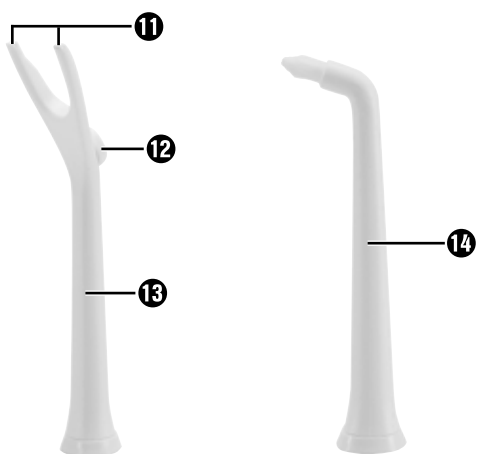**A****B**

# Contents

<b>Introduction</b> .....	<b>2</b>
<b>Intended use</b> .....	<b>2</b>
<b>Package contents</b> .....	<b>2</b>
<b>Appliance description</b> .....	<b>2</b>
<b>Technical details</b> .....	<b>2</b>
<b>Safety information</b> .....	<b>3</b>
<b>Setting up</b> .....	<b>4</b>
<b>First use</b> .....	<b>4</b>
<b>Fitting/removing the brush head/interdental brush/ floss attachment</b> .....	<b>4</b>
<b>Operation</b> .....	<b>5</b>
<b>Cleaning your teeth with the brush head</b> .....	<b>5</b>
<b>Cleaning with the interdental brush</b> .....	<b>6</b>
<b>Cleaning with the floss attachment</b> .....	<b>7</b>
<b>Cleaning and care</b> .....	<b>8</b>
<b>Travel lock</b> .....	<b>9</b>
<b>Disposal</b> .....	<b>9</b>
<b>Ordering replacement parts</b> .....	<b>10</b>
<b>Kompernass Handels GmbH warranty</b> .....	<b>11</b>
Service .....	12
Importer .....	12

## Introduction

Congratulations on the purchase of your new device.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the range of applications specified. Please also pass these operating instructions on to any future owner.

## Intended use

This appliance is intended exclusively for personal dental care at home. This appliance is intended solely for use in private households. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. Do not use it for commercial purposes.

## Package contents

Handset

8 brush heads

Charger

Operating instructions

### NOTE

Check contents for completeness and damage immediately after unpacking. Contact customer service if necessary.

Remove all packaging materials and plastic films from the appliance.

## Appliance description

Figure A:

- 1 Brush head
- 2 Drive shaft
- 3 Lid
- 4 Brush holders
- 5 Charger
- 6 Charge indicator
- 7 "Pulse" LED (yellow)
- 8 "Soft" LED (green)
- 9 "Clean" LED (orange)
- 10 On/Off switch

Figure B:

- 11 Guide grooves
- 12 Fixture
- 13 Floss attachment
- 14 Interdental brush

## Technical details

Power supply: 100–240 V ~ (alternating current), 50/60 Hz

Power consumption: 2 W

Rechargeable battery: Li-ion battery, 600 mAh, 3.7 V

Protection type (charger): IP X7 (protected against the ingress of water during temporary immersion)

Protection class: II / 

Charging time: approx. 14 - 16 hours (first charge)

Protection type (handset): IP X7 (protected against the ingress of water during temporary immersion)

## Safety information

### **⚠ DANGER – RISK OF ELECTRICAL SHOCK!**

- ▶ To avoid damage, do not allow liquids to enter the appliance and do not submerge it in water.
- ▶ Ensure that no foreign objects enter the appliance.
- ▶ In the event of malfunctions or obvious damage, please contact our customer service department.
- ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way. Danger of electric shock!
- ▶ The mains cable of this appliance cannot be replaced. If the cable is damaged, the whole appliance must be scrapped.
- ▶ Never touch the charger if it has fallen into water. Remove the plug from the wall socket first! Danger of electric shock!

 Do not submerge the handset or the charger in water or other liquids.

### **⚠ WARNING – RISK OF INJURY!**

- ▶ Toothbrushes may be used by children and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge, if they are supervised or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance should not be carried out by children, unless they are under supervision.
- ▶ Do not allow children to play with the appliance.

### **CAUTION – PROPERTY DAMAGE!**

- ▶ The electric toothbrush may only be recharged with the supplied charger. Otherwise, the appliance could be irreparably damaged!
- ▶ Operate the appliance only with the original accessories supplied.

## Setting up

### NOTE

- ▶ In damp rooms, especially in the vicinity of water (bath), the charger ⑤ must be set up so that it cannot accidentally fall into the water.

- 1) Place the charger ⑤ on a clean and level surface.  
Ensure that the charger ⑤ cannot fall into water, such as into a bath or sink.

## First use

### NOTE

Before starting to use the appliance, check that it is in fault-free condition. You must not use the appliance if it has any faults.

- 1) Insert the power plug into the socket.
- 2) Before the first use, clean the new brush head ① / the interdental brush ⑭ / the floss attachment ⑬ with clear water.
- 3) Open the lid ③ and place the replacement brush heads, the interdental brush ⑭ and/or the floss attachment ⑬ into the brush holders ④ as required. Close the lid ③.
- 4) Place the handset on the charger ⑤ so that it slides onto the pin on the base.
- 5) Allow the handset to charge for around 14–16 hours before using it for the first time. The charge indicator ⑥ on the front of the handset lights up red during the charging process.
- 6) Use the handset repeatedly without charging the battery until the performance becomes noticeably weaker or the toothbrush stops working and the charge indicator ⑥ flashes red. Then recharge the battery for 14–16 hours. When the appliance is fully charged, the charge indicator ⑥ lights up green.
- 7) Charge and discharge the toothbrush completely 3 times. At that point, the battery will have reached its full capacity.
- 8) You can now leave the handset on the charger ⑤ after every use, so as to always have full power available.

## Fitting/removing the brush head/interdental brush/floss attachment

- To remove the respective attachment pull it upwards and off the drive shaft ②.
- When fitting the attachment, ensure that the drive shaft ② on the handset slides precisely into the socket on the attachment.  
Push the attachment down until it completely covers the drive shaft ②.

### NOTE

We recommend only using original attachments. You can order replacement attachments at any time (see section "Ordering replacement parts").

- ▶ You can also use the replacement brush heads of the model Oral B® Sonic Complete\* for this appliance.

\*Oral-B is a registered trademark of Procter & Gamble.

## Operation

This toothbrush has 3 different power levels:

- “Clean”: normal speed, orange “Clean” LED ⑨  
For cleaning the teeth.
- “Soft”: slow speed, green “Soft” LED ⑧  
For cleaning the teeth of people with sensitive gums.
- “Pulse”: slow speed with short breaks, yellow “Pulse” LED ⑦  
For massaging the gums.

- Press the on/off switch ⑩ once to activate the “Clean” mode. The “Clean” LED ⑨ lights up.
- Press once more to switch to “Soft” mode. The “Soft” LED ⑧ lights up.
- Press it a third time to switch to “Pulse” mode. The “Pulse” LED ⑦ lights up.
- To switch the toothbrush off, press the on/off switch ⑩ once more.

### NOTE

If you operate one of the modes for longer than 10 seconds and then press the On/Off switch ⑩, the toothbrush will switch off.

Next time you switch it on, the toothbrush will automatically start in the mode you were using last.

To change the mode again, press the On/Off switch ⑩ repeatedly until the appliance switches to the required mode as described above.

## Cleaning your teeth with the brush head

### ⚠ WARNING!

The unfamiliar stimulation may cause light gum bleeding at first. This is normal. If this continues for more than two weeks, you should consult your dentist.

- 1) Place the brush head ① on the handset.
- 2) Moisten the brush under flowing water.
- 3) Put a little toothpaste onto the wet brush.

### NOTE


To avoid being sprayed, do not switch on the appliance until it is in your mouth.

- 4) Use the on/off switch ⑩ to select the desired operating mode.

### NOTE

This toothbrush has a 30-second interval timer, so that all four jaw quadrants can be cleaned evenly.

- 5) Divide your mouth into four parts (quadrants): upper left, upper right, lower left and lower right.

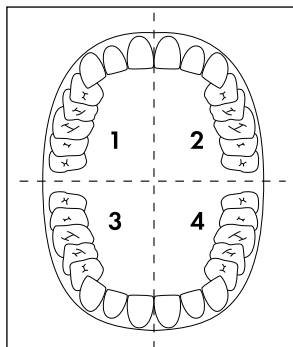- 6) Start in one quadrant:  
Start with the brush head **1** on your gums, and move the brush upwards/  
downwards from the gums to the crown of the tooth.  
Continue in this manner on both the inside and the outside of the tooth.  
Start again at the gums for each individual tooth.  
Do not forget to brush the back of the rear molars.  
To clean the chewing surfaces, guide the brush over them with a little  
pressure.
- 7) After 30 seconds, the toothbrush stops briefly to signal that you should  
start cleaning the next quadrant. Clean each jaw quadrant using the same  
technique.
- 8) After two minutes, the toothbrush stops briefly 2 times. This informs you  
that you have complied with the dental recommendation of two minutes  
cleaning.
- 9) Press the on/off switch **10** repeatedly until the toothbrush stops.
- 10) Then take the toothbrush out of your mouth.
- 11) After use, rinse the brush head **1** thoroughly under running water.
- 12) Remove the brush head **1** from the handset.  
To clean the drive shaft **2** and the handset, hold them under running wa-  
ter for a short while. Hold the handset so that the drive shaft **2** is pointing  
downwards at an angle.

**NOTE**

Massage your gums regularly. This promotes blood circulation and helps to keep them healthy. Set the mode to "Pulse" and guide the moistened brush over your gums with a light pressure.

## Cleaning with the interdental brush

You can use the interdental brush **14** to clean the interdental spaces, difficult-to-access parts on the teeth, or the chewing surfaces of your molars.

- 1) Attach the interdental brush **14** to the hand element.
- 2) Moisten the interdental brush **14** under running water.
- 3) Put a little toothpaste onto the wet interdental brush **14**.

**NOTE**

To avoid splashes, do not switch on the appliance until it is in your mouth.

- 4) Use the On/Off switch **10** to select the desired operating mode and clean the required area.
- 5) Press the On/Off switch **10** repeatedly until the appliance stops.
- 6) Then take the interdental brush **14** out of your mouth.
- 7) After use, clean the interdental brush **14** thoroughly under running water.
- 8) Remove the interdental brush **14** from the hand element.  
To clean the drive shaft **2** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **2** is pointing downwards at an angle.

**NOTE**

The interdental brush **14** is ideal for cleaning braces or false teeth.

## Cleaning with the floss attachment

You can use the floss attachment **13** to clean interdental spaces with the aid of dental floss.

- 1) Place the floss attachment **13** on the hand element.
- 2) Take around 15 cm of standard commercial dental floss and wrap one end around the fixture **12** 2-3 times.

### NOTE

If the dental floss slips off the fixture **12**, you can tighten the screw with a Phillips screwdriver.

- 3) Pass the dental floss tightly through the guide grooves **11** on the floss attachment **13**.
- 4) Then wrap the other end around the fixture **12** 2-3 times so that the dental floss is stretched taut and does not come undone while cleaning the interdental spaces.
- 5) Guide the dental floss between two teeth with the toothbrush switched off.
- 6) Switch the appliance on. When cleaning with the floss attachment **13**, use the settings "Clean" or "Soft".
- 7) Clean the interdental space with careful up-and-down movements.
- 8) Clean all the interdental spaces as described above.
- 9) When you are finished, switch off the toothbrush and take it out of your mouth.
- 10) Remove the dental floss from the fixture **12** and dispose of it.
- 11) After use, clean the floss attachment **13** thoroughly under running water.
- 12) Remove the floss attachment **13** from the hand element.  
To clean the drive shaft **2** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **2** is pointing downwards at an angle.

## Cleaning and care

### **DANGER – RISK OF ELECTRICAL SHOCK!**

- ▶ Always remove the plug of the charger **5** from the power socket before cleaning the appliance! There is a risk of electric shock!

 Do not immerse the handset or the charger **5** in water or other liquids!

### **WARNING – PROPERTY DAMAGE!**

- ▶ Do not use detergents or solvents. They can damage the plastic surfaces.

### **NOTE**

Never open the housing of the charger **5** or the toothbrush! They cannot be resealed in a watertight manner.

- 1) Dismantle the charger **5** for cleaning:
  - Hold the lower part of the charger **5** firmly and tilt the top part, which contains the brush holders **4**, towards the rear.
- 2) Clean the upper part of the charger **5** under running water. If any toothpaste residue in the brush holders **4** cannot be removed, soak the upper part of the charger **5** in warm water and then rinse out the brush holders **4**.
- 3) Clean the lower part of the charger **5** only with a soft, slightly damp cloth.
- 4) Remove the brush head **1** from the handset. After use, clean the brush heads **1** thoroughly under running water.  
We recommend cleaning the brush heads **1** in a standard commercial denture cleaning fluid every 2–4 weeks.  
The brush heads **1** are equipped with so-called indicator bristles (the outer, light-blue bristles). These bristles lose their colour to show that the brush head **1** needs to be changed. Change the brush head **1** when the indicator bristles lose their colour.
- 5) Clean the interdental brush **14** after each use under running water. You should replace the interdental brush **14** after about 3 months of regular use, or when the bristles start to bend outwards.
- 6) Clean the floss attachment **13** under running water after every use. If it is damaged, however, it must be replaced.
- 7) To clean the drive shaft **2** and the handset, hold them under running water for a short while. Hold the handset so that the drive shaft **2** is pointing downwards at an angle.

## Travel lock

This toothbrush is equipped with a travel lock. You can lock the on/off switch **10** to prevent accidental starting of the toothbrush:

Press and hold the on/off button **10** until the charge indicator **6** flashes green 5 times.

The on/off switch **10** is now locked.

To unlock the switch, press and hold the on/off button **10** until the charge indicator **6** lights up once for a few seconds.

### NOTE

If you try to switch on the toothbrush while the travel lock is activated, the charge indicator **6** will flash 5 times. The toothbrush will not start.

## Disposal

**Never dispose of the appliance in your normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.**

Dispose of the device via an approved disposal company or your municipal waste facility. Please comply with all applicable regulations. Please contact your waste disposal centre if you are in any doubt.


Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

### NOTE

- ▶ The integrated rechargeable battery in this appliance cannot be removed for disposal.
- ▶ To prevent risks, the removal or replacement of the battery may only be carried out by the manufacturer or the manufacturer's customer service department or a similarly qualified person.
- ▶ When disposing of the device, indicate that the device contains batteries.

The packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

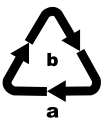Dispose of the packaging in an environmentally friendly manner.

Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

- 1-7: Plastics,
- 20-22: Paper and cardboard,
- 80-98: Composites

## Ordering replacement parts

You can order the following replacement parts for the product NSZB 3.7 B3:

- ▶ 8-pack of toothbrush heads

- ▶ 2-piece set  
(interdental brush and floss attachment)

Order the replacement parts via the service hotline (see section "Service") or simply visit our website at [www.kompernass.com](http://www.kompernass.com).

### NOTE

- ▶ Have the IAN number of the appliance ready (can be found on the cover of these operating instructions) when you place your order.

## Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

### Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

### Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

### Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

## Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on [www.lidl-service.com](http://www.lidl-service.com).

## Service

### **GB** Service Great Britain

Tel.: 0871 5000 720 (£ 0.10/Min.)

E-Mail: [kompernass@lidl.co.uk](mailto:kompernass@lidl.co.uk)

### **IE** Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: [kompernass@lidl.ie](mailto:kompernass@lidl.ie)

IAN 303047

## Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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